



Pacifiers

Medical Device - Hospital Use



LIFE STARTS WITH US

Babies are extremely sensitive patients whose nervous, immune and digestive systems are still maturing.

Non-nutritive sucking is essential for the proper development of the premature baby and for the management of newborns' pain. The priority is to support healthcare professionals and parents in providing care to these small sensitive patients in a hygienic and completely safe manner.



Support the **development of oral skills** in premature babies and **manage pain** in newborns through non-nutritive sucking

Beldico single-use pacifiers are medical devices ready for use and designed to **administer a sugary analgesic solution** or **satisfy the non-nutritive sucking needs** of premature, newborns and infants in a hospital setting.

The shape and material of Beldico pacifiers have been carefully selected for **the comfort and safety** of premature babies, newborns and infants as well as for convenience for caregivers and parents.

Beldico pacifiers, sterile or bacteriologically clean, individually packed, meet the needs of **neonatology, maternity** and **paediatric** departments.

0%

BPA/BPS-free

0%

Latex-free

0%

Phthalates-free



Food compatible




Tasteless and odourless

Considered to be an essential care for premature babies' development, **non-nutritive suction significantly decreases the length of hospital stay in preterm infants^[1].**

^[1] Pinelli, J., & Symington, A. (2000). How rewarding can a pacifier be? A systematic review of nonnutritive sucking in preterm infants.

Designed for the needs of preterm and newborn babies

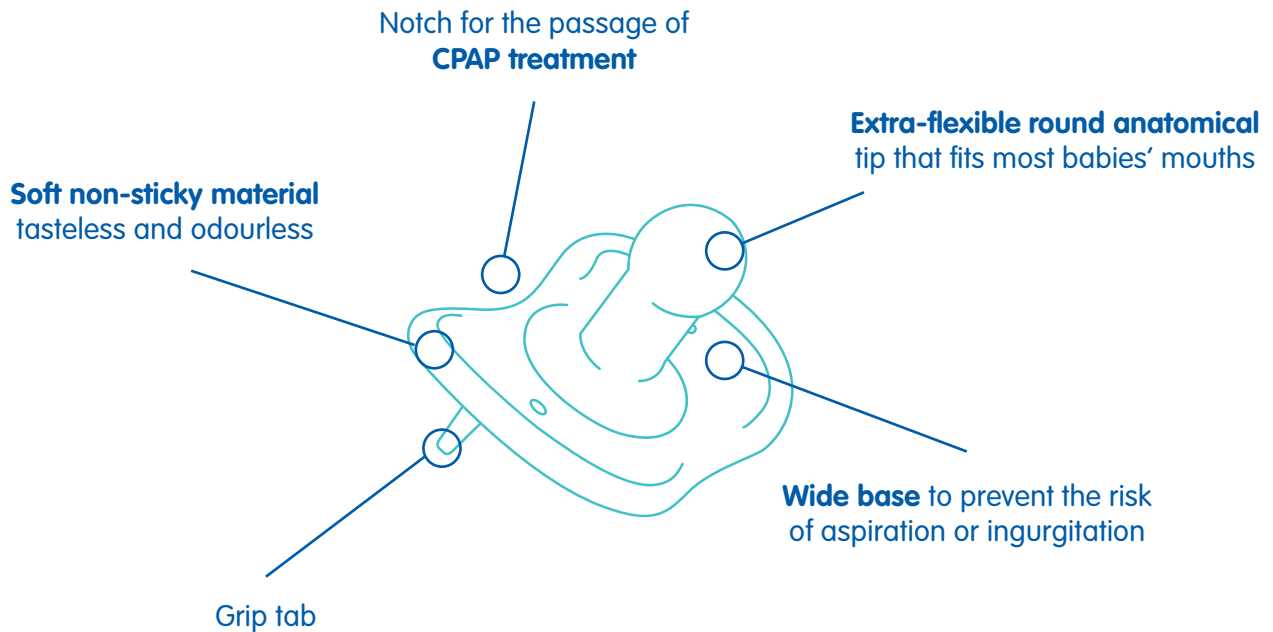

Beldico pacifiers are used in the following cases

- **ASSOCIATED WITH A CPAP TREATMENT OR NASAL CPAP TREATMENT:** for premature babies and newborns.
- **ASSOCIATED WITH ENTERAL FEEDING:** during enteral nutrition to promote the association between sucking and nutrition.
- **FEEDING TRANSITION SUPPORT:** to facilitate the transition from enteral nutrition to bottle feeding or breast-feeding after stopping enteral nutrition.
- **SUPPORT FOR THE WEANING PROCESS:** after the disconnection from intubation or, other respiratory treatments.
- **SUPPORT FOR PAIN MANAGEMENT:** during painful procedures in premature and term newborns, the pacifier is a simple, non-invasive, effective, adjuvant in the management of pain and stress. It can be immersed in 30% glucose or sucrose solutions before being offered to the baby.

2 sizes available

Adapted according to the baby's weight



Premature baby
1 to 2 kg



Newborn
2 to 6 kg

Suited to hospital practices

- **SAFE:** Beldico pacifierstare **made in a single piece**, their shape and their size are tailored to the baby's weight to protect against the risks of aspiration or of ingurgitation. Made from food-grade materials, Beldico pacifiers are certified latex-free, phthalate-free, and free of BPA and BPS.
- **HYGIENIC:** Beldico pacifiers are **sterile or bacteriologically clean** and packed individually, guaranteeing their hygiene until use by the infant.
- **ADAPTED:** Thanks to their **extra-flexible, anatomical round tip**, Beldico pacifiers fit the majority of babies' mouths. Their soft, tasteless and odourless material is appreciated by babies.

Standard features

Name	Single-use pacifiers
Material	Pacifier: Thermoplastic elastomer (TPE) Guaranteed latex- free, phthalate- free, BPA- free, BPS- free
Type of treatment	Sterilization SAL 10 ⁻⁶ by ETO or disinfection SAL 10 ⁻² by gamma radiation
European Directive 93/42/EEC	Medical Device class I (bacteriologically clean) Medical Device class Is (sterile)
Standards	EN 1400 Child use and care articles. Soothers for babies and young children. Safety requirements and test methods (section 10 only) ISO 14644-1 & 2 - clean rooms and associated controlled environment EN ISO 10993-1 Biological evaluation of medical devices ETO: EN ISO 11135 Sterilization of health-care products - Ethylene oxide GAMMA: EN ISO 11137-1 & 2 - sterilization of health care products - radiation

Beldico develops, produces and markets medical devices for feeding babies and focuses particularly on the specific needs of newborns and premature babies.

Bottles - breastshields - teats - pacifiers
Breast milk expression and handling
Breast milk or formula warming

Manufacturer

BELDICO
Rue André Feher 3
BE-6900
Marche-en-Famenne
BELGIUM
Tel : 0032 (0)84 31 37 80
www.beldico.be

Distributor