


Order information

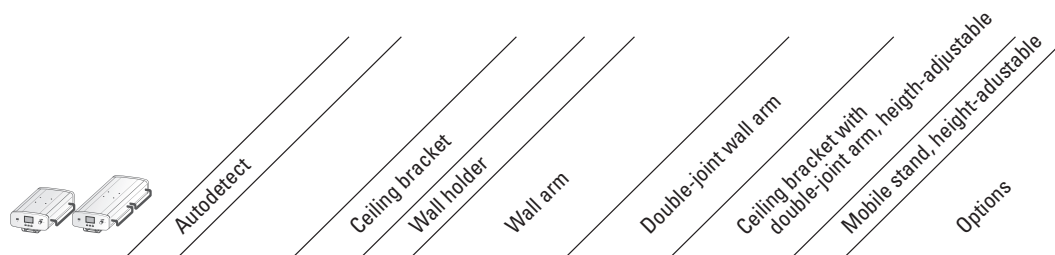

Order No.	CERAMOTHERM 3100 4,8 kg	CERAMOTHERM 3200 6,8 kg	Automatic distance detection	Weight (kg) complete	Torque (Nm)	Static, distance between patient pad and solid ceiling max. 2600 mm	Static	320 mm length	480 mm length	600 mm length	R1 = 480 mm, R2 = 480 mm	Height-adjustable, R1 = 750 mm, R2 = 800 mm	Distance between patient pad and solid ceiling max. 2600 mm	Distance between patient pad and solid ceiling above 2600 mm	Vertical adjustment ± 45° WY1822	Vertical extension by 250 mm WY1846	Wall arm bracket support M0369	Extension arm R = 950 mm Wall: WY1847, Ceiling: WY1848
WY3101	•			8,8	27			•							•	•	•	
WY3102	•			9,2	39				•						•	•	•	
WY3103	•			9,4	48					•					•	•	•	
WY3104	•			12,8	90						•						•	
WY3107	•		•	17,8	181							•					•	•
WY3108*	•		•	27,4	171								•					
WY3112	•			22,3	–									•				
WY3115	•			7	17			•										
WY3116	•			6,3	11			•L										
WY3117*	•			12,3	13,3	•												
WY3201		•		10,9	35			•L								•	•	•
WY3202		•		11,2	50				•							•	•	•
WY3203		•		11,4	62					•						•	•	•
WY3204		•		14,8	111						•						•	•
WY3207		•		19,8	214							•					•	•
WY3208*		•		29,9	202								•					•
WY3212		•		23,3	–									•				
WY3216		•		8,3	16			•L										
WY3217*		•		14,3	59	•												
B0329*				30 _{kg/m}										•				
M0369				2,1														

L Lateral mounting

* **When ordering please indicate dimensions:** Floor to top line of patient pad
 Floor to bottom line of suspended ceiling
 Floor to bottom line of solid ceiling

Mounting information

Wall mount with flexible arrangement to the patient pad

Mounting in solid brickwork, lime sandstone, concrete, according to the national construction standards.

In agreement with the architect or engineering office cavity and light-weight construction walls must be reinforced as follows: Height 600, width according the distance between the posts. Material: Sheet steel min. 2 mm thick or multiplex plate min. 35 mm thick, screwed in place from post to post.

As an alternative, for double-planked plasterboard walls (2 x 12.5 mm) the external wall arm bracket support order No. M0369 may be used. In that case expansion anchors Hilti HDD-S 8 or equivalent must be used and the wall arm bracket must be screwed in place on the wall arm bracket support by min. 2 x M6 (8.8) screws.

Order No. WY3101, WY3102/3202, WY3103/3202, WY3104/3204, WY3107/3207

Ceiling mounts

Mounting at the solid ceiling

Static: Distance between patient pad and bottom line of solid ceiling max. 2600 mm respectively distance between solid ceiling and suspended ceiling max. 500 mm.

Order No. WY3117, WY3217

Height-adjustable: Distance between patient pad and bottom edge of solid ceiling max. 2600 mm respectively distance between solid ceiling and suspended ceiling max. 500 mm.

Order No. WY3108, WY3208

In case of distance between patient pad and bottom line of solid ceiling above 2600 mm respectively distance between solid ceiling and suspended ceiling above 500 mm a connection flange B0329 is required in addition.